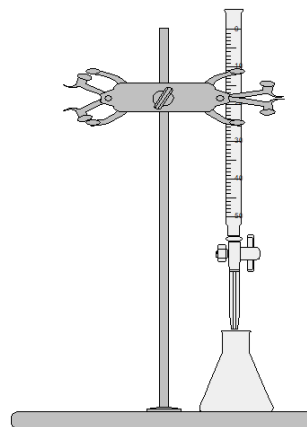

Experiment 14 – Acid/Base Titrations

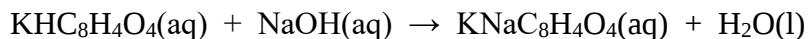
Discussion

This experiment demonstrates an analytical technique known as titration, where a solution is delivered from a buret until it completely consumes another solution in a flask. Consider the following:

Acid-base titrations are an example of **volumetric analysis**, a technique in which one solution is used to analyze another. The solution used to carry out the analysis is called the **titrant** and is delivered from a device called a **buret**, which measures the volume accurately. The point in the titration at which enough titrant has been added to react exactly with the substance being determined is called the **equivalence point** (or **stoichiometric point**). This point is often marked by the change in color of a chemical called an **indicator**. The titration set-up is illustrated in the schematic shown.



First, the concentration of a base solution (standard) will be determined; this data will be used to determine the concentration of unknown acid solutions. To standardize the NaOH solution, it will be reacted with potassium hydrogen phthalate, $\text{KHC}_8\text{H}_4\text{O}_4$ (abbreviated KHP). The molar mass of KHP is 204.2 g/mol. The reaction of KHP with NaOH is known to be:



Starting with a known mass of KHP and recording the volume of NaOH required to reach the endpoint, the molarity of the base can be determined. The indicator to be used, phenolphthalein, is colorless in acidic solution and rosy pink when slightly basic. Even though the endpoint is slightly overshoot in order to make the color change, the goal is to use as little excess base as possible. Therefore, the titrated solution should be very pale pink, not bright rosy red, at the endpoint.

Once the concentration of the NaOH solution is known, one could determine the concentration of an acid solution by doing an additional titration.

Procedure

1. Measure out 1.000 to 1.200 g of KHP and add to a 125 mL Erlenmeyer flask.
2. Add approximately 30 mL of D.I. water to the flask. If some KHP is sticking to the walls of the flask, rinse it down with D.I. water from a wash bottle.
3. Take a clean, dry 100 mL beaker to the hood and obtain approximately 60 mL of NaOH base. Be sure to keep the beaker covered with a watch glass when not in use.
4. Obtain a clean 50 mL buret. Carefully fill the buret with base, making sure that no air bubbles are present. Run some of the base solution through the buret tip to remove the air pocket in the tip.
5. Record the initial buret reading in the data section. When you read a buret, the line of sight must be level with the BOTTOM of the meniscus to avoid error. The top of the buret reads 0.00 mL; the bottom reads 50.00 mL. Notice the numbers increase going down. Take note of this when reading the numbers. Your instructor will demonstrate.
6. Add 2 to 3 drops of phenolphthalein indicator solution to the 125 mL Erlenmeyer flask containing KHP and water. Swirl your acid solution until the KHP is completely dissolved before moving on to the next step of the procedure.
7. Place the flask under the tip of the buret. A piece of white paper under the flask makes it easier to see the pale pink color at the endpoint. Open the valve and allow base to drip from the buret into the flask. Swirl continually to mix the solutions. As you get closer to the endpoint, the solution will begin to show pink color that goes away when you mix. Slow the rate of base addition to one drop at a time, mixing the solutions well after every drop. If you splash the solution up onto the sidewalls of the flask, spray a stream of water from your wash bottle over the inside of the flask. The extra water that mixes into your acid sample will not affect your results. Once the addition of ONE drop of base changes the solution from colorless to pale pink, close the buret valve and make certain that the pale pink color lasts for at least 30 seconds while swirling. If so, record the final buret reading. If not, carefully add one more drop of base from the buret valve until the pale pink color persists for 30 seconds of swirling or longer. If at the end of your trial the color is bright rosy red, you have overshoot the endpoint. Make a note in your data if you overshoot.
8. Discard the titrated solution into the sink.
9. You may be instructed to re-fill the buret to near the 0.00 mL line with NaOH solution and stopper to save the solution for the next lab section.
10. Calculate the molarity of the base solution from your trial. Calculate the average molarity of the base from the class data.

Name: _____

Section: _____

Data and Calculations for Experiment 14

Be sure to record all numerical values to the correct significant figures and include units!

Mass of flask and KHP	
Mass of empty flask	
Mass of KHP	
Initial buret reading	
Final buret reading	
Volume of base used	

CALCULATIONS – Show your work!

1. Moles of acid (KHP, Molar mass = 204.2)

2. Moles of base used to neutralize acid

3. Molarity of base (NaOH)

4. *Class Average* Molarity of Base: